

## Instrument Installation and Operational Services

Installation Qualification (IQ) and Operational Qualification (OQ) services are performed by our Field Service Engineers and are critical to ensuring that your instrument is properly installed, and operating according to specifications. Luminex offers standard IQ/OQ for new instruments, and Yearly OQ to confirm system reliability of any xMAP® instrument in your laboratory.

**Each service includes a final protocol document** which can be used to fulfill a broad range of regulatory and quality requirements.

### Why are these services needed?

- You will benefit from the expertise of engineers that know your instrument inside and out.
- IQ/OQ protocols are the evidence you need to prove that your instrument is properly installed, and operating according to specifications.
- Yearly OQ provides additional assurance that the system will produce precise and reliable assay results, even after being repaired or moved to a different location in the lab.

### IQ/OQ

The IQ service confirms that the instrument, operating software, and other accessories have been supplied as ordered and that the system has been installed per Luminex guidelines.


#### During the IQ, we will:

- Verify that the installation site meets the environmental requirements for the instrument (bench space, power, temperature, etc.)
- Review the hazards and safety precautions pertinent to the use of the Luminex system
- Provide installation checklists for the software, and hardware (including tubing, cables, and accessories)
- Document preventative maintenance (PM) requirements for inclusion in the PM schedule

Once the system is installed, the OQ service will be conducted to demonstrate that the system is functioning according to its operational specifications.

#### The OQ will:

- Test all components of the Luminex system for proper functionality
- Challenge instrument specifications using calibration and performance verification operations

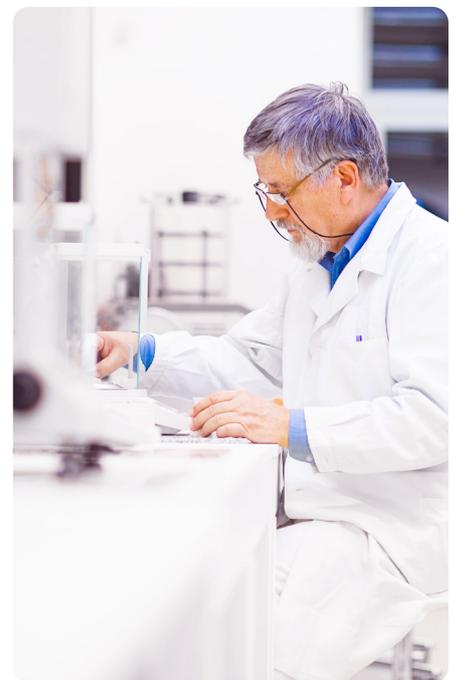

## Yearly OQ

The Yearly OQ service verifies lab environmental fitness for optimal system performance and is highly recommended to confirm system reliability on an annual basis or after a renovation event or instrument repair. Yearly OQ offers additional services which are not a part of routine PM.

### Yearly OQ includes:

- All standard OQ verifications, plus some additional IQ parameters
- Identifying environmental factors that may affect system performance such as—electrical changes in the lab or other parts of the facility, HVAC maintenance, facility room renovations, or changes to bench top and capital equipment configurations within the lab
- Checks to ensure that maintenance has been adequately performed by lab personnel

When combined with routine PM, Yearly OQ can help to ensure that your system will produce precise and reliable assay results.

## Ordering Information

System	Initial IQ/OQ Part Number	Yearly OQ Part Number
MAGPIX®	VP-MAGPIX-XPON-IQOQ	VP-MAGPIX-XPON-OQ
Luminex® 200™	VP-LX-XPON-IQOQ	VP-LX-XPON-OQ
FLEXMAP 3D®	VP-FM3D-XPON-IQOQ	VP-FM3D-XPON-OQ
xMAP INTELLIFLEX®*	VP-INTELLIFLEX-IQOQ	VP-INTELLIFLEX-OQ

All IQ/OQ and Yearly OQ services are performed solely by experienced Luminex Field Service Engineers. Please contact your Luminex representative for information on how to order this service.

**Luminex**  
A DiaSorin Company

[luminexcorp.com/eu/installation-qualification-iq-operational-qualification-oq-services/](https://www.luminexcorp.com/eu/installation-qualification-iq-operational-qualification-oq-services/)

\*For Research Use Only. Not for use in diagnostic procedures.

©2012-2023 Luminex Corporation. A DiaSorin Company. All rights reserved. Luminex, xMAP, MAGPIX, and FLEXMAP 3D are trademarks of Luminex Corporation, registered in the U.S. and other countries. 200 is a trademark of Luminex Corporation.

[www.luminexcorp.com](https://www.luminexcorp.com)

HEADQUARTERS

UNITED STATES

+1.512.219.8020

info@luminexcorp.com

EUROPE

+31.73.800.1900

europe@luminexcorp.com

CANADA

+1.416.593.4323

info@luminexcorp.com

CHINA

+86.21.8036.9888

infocn@luminexcorp.com

JAPAN

+81.3.5545.7440

infojp@luminexcorp.com

FL517532